


These rapid drying moisture extraction ovens have variable forced air extraction.

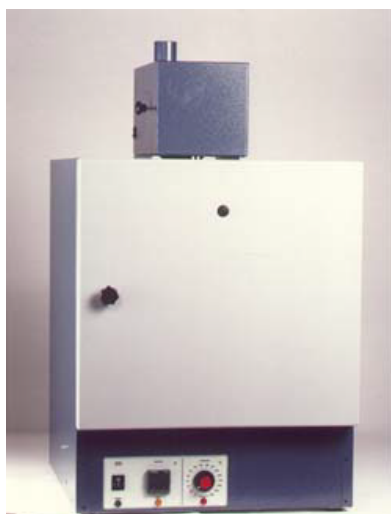

The exterior is constructed from sheet steel finished in easy clean stove powder paint. The interior is stainless steel (SS) and the shelves are chrome plated steel wire grids supported on fixed shelf runners.

The door and chamber are fitted with silicone seals for excellent insulation from heat and are easy to clean.

Incoloy sheathed heating elements are controlled by either a dial controller or digitally.

There are 10 models in the range with capacities from 50 to 350 litres.

Standard Specifications:

- Temperature range 40C to 250C
- Stainless steel interior
- Fan assisted circulation
- Output vent for 50mm diameter tubing
- Direct reading hydraulic temperature control
- Safety thermostat with warning indicator

Options Include:

- Extract boost fan & variable output vent
- Digital control and readout (DIG)
- Dial thermometer readout
- Programmable timers
- Adjustable shelf ladders

Model	Litres	Shelf Positions	No of Shelves	External mm H x W x D	Internal mm H x W x D	Max Load Height mm	Power watts
GPME50	50	3	2	615 x 590 x 422	325 x 490 x 315	240	1000
GPME75	75	3	2	623 x 590 x 567	333 x 490 x 460	248	1250
GPME100	100	4	3	735 x 590 x 567	445 x 490 x 460	360	1250
GPME120	120	4	3	735 x 590 x 659	445 x 490 x 552	360	1500
GPME150	150	5	3	845 x 590 x 659	555 x 490 x 552	470	1500
GPME180	180	6	4	957 x 590 x 659	667 x 490 x 552	582	1750
GPME210	210	6	4	957 x 590 x 749	667 x 490 x 642	582	2000
GPME240	240	8	4	1178 x 590 x 659	888 x 490 x 552	803	2250
GPME280	280	6	4	941 x 750 x 769	651 x 650 x 662	566	2250
GPME350	350	7	4	1104 x 750 x 769	814 x 650 x 662	729	2500

Standard voltage 230v, 50/60Hz Other voltages to order